

ALOVA[®] CONTOURED CUSHION



VISCO ELASTIC MOULDED CUSHION WITH A COCCYX CUT-OUT. FOR PERSONS UP TO 200 KG (31 STONE).

- **Efficiency** : the ALOVA[®] visco elastic material increases the contact area to better distribute and reduce maximum pressure values.
- **Safety** : reduced risk of slipping forward thanks to the tilted seating position. Integral cover with anti slip base.
- **Comfort** : thermoregulation thanks to Promust PU properties.
- **New** : the coccyx cut-out reduces the risk of pressure ulcers on the coccyx and ischiatic regions (cushion more than 6 cm thick under ischions).



ALOVA® CONTOURED


Efficiency :

- Coccyx cut-out: improved prevention of pressure ulcers in sacrum and coccyx regions.
- Minimum thickness of the cushion: 6 cm
Improved prevention of pressure ulcers in ischiatic regions.

Ergonomics :

- 16 sizes and 3 heights: maximum patient weight up to 200 Kg (31 stone).
- Pelvic support to reduce slipping forward.

- Anti slip base.
- Removable cover.
- Handle on the left base of the cushion: less risks of tearing and more hygienic.

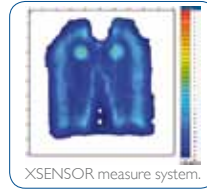

Evaluation using pressure XSENSOR for a patient of 80 Kg and 172 cm.

ALOVA® ERGONOMIC CUSHION

- The pads follow the movements of the body and ensure air circulation.
- One single size: 43x43x9 cm.
- Thermoregulation thanks to the vapour permeable cover and excellent air circulation between the pads.
- Immersion into Alova foam to reduce pressure points.



TECHNICAL SPECIFICATIONS

DESIGNATION	REFERENCES	WIDTH CM	DEPTH CM	HEIGHT CM MINI / MAXI*	PATIENT WEIGHT KG (STONE)	COVER	FIRE STANDARDS
ALOVA® Ergonomic Cushion **	VCAE/M	43	43	9	25 to 100 Kg (3 to 15 st)	Integral Promust PU cover: Jersey coated PU. Anti slip base (PU/PVC material)	EN 597-1 EN 597-2
ALOVA® Contoured Cushion Standard size with 2 covers **	VCAG4343	43	43	6 / 9	30 to 100 Kg (4 to 15 st)	Integral cover in Promust PU (Jersey coating polyurethane) black on the top and silver on the sides.	
ALOVA® Contoured Cushion multi-sizes with 2 covers **	VCAG4340	43	40	6 / 9	30 to 90 Kg (4 to 14 st)		
	VCAG4540	45	40	6 / 9	30 to 120 Kg (4 to 18 st)		
	VCAG3841	38	41	6 / 9	30 to 90 Kg (4 to 14 st)		
	VCAG4041	40	41	6 / 9	30 to 90 Kg (4 to 14 st)		
	VCAG4341	43	41	6 / 9	30 to 100 Kg (4 to 15 st)		
	VCAG4043	40	43	6 / 9	30 to 90 Kg (4 to 14 st)		
	VCAG4544	45	44	6 / 9	30 to 120 Kg (4 to 18 st)		
ALOVA® Contoured Cushion multi-sizes XL with 2 covers **	VCAG4343XL-E	43	43	9 / 12	30 to 110 (4 to 16 st)	Anti slip base PU/ PVC equipped with an handle.	
	VGAG5044XL-E	50	44	9 / 12	30 to 135 Kg (4 to 21 st)		
ALOVA® Contoured Cushion multi-sizes XXL with 2 covers **	VCAG5250	52	50	7 / 10	120 to 160 (18 to 25 st)		
	VCAG5450	54	50	7 / 10	120 to 160 (18 to 25 st)		
	VCAG5650	56	50	7 / 10	120 to 200 Kg (18 to 31 st)		
	VCAG6050	60	50	7 / 10	120 to 200 Kg (18 to 31 st)		
ALOVA® Contoured with Pelvic support with 2 covers	VCAG3841BP	38	41	6 / 9	30 to 90 Kg (4 to 14 st)		
	VCAG4041BP	40	41	6 / 9	30 to 90 Kg (4 to 14 st)		
	VCAG4341BP	43	41	6 / 9	30 to 100 Kg (4 to 15 st)		
	VCAG4343BP	43	43	6 / 9	30 to 100 Kg (4 to 15 st)		
	VCAG4544BP	45	44	6 / 9	30 to 120 Kg (4 to 18 st)		
	VCAG5044BP	50	44	6 / 9	30 to 135 Kg (4 to 21 st)		
	VCAG4347BP	43	47	6 / 9	30 to 120 Kg (4 to 18 st)		

* minimum height = minimum thickness of the cushion under the ischions / maximum height = total height of the cushion.

** If you want to order a cushion with only one cover, please add **-E** at the end of the reference

Cover compatible with detergent and disinfecting products for hospital use.

Care of Promust PU Cover



Cushion guarantee: 3 years



WINCARE
g r o u p

by

ASKLÉSANTÉ

ASKLÉ SANTÉ Groupe WinnCare – 200 rue Charles Tellier – ZI de Grézan – 30034 Nîmes Cedex 1 France

Tél : +33 (0)4 66 02 15 15 – Fax : +33 (0)4 66 02 15 28 – E-mail : contact@askle.com

www.winnCare.fr

ASKLÉ SANTE reserves the right to modify technical characteristics without prior notice - non contractual photos.